



The Worlds Largest manufacturer of x-ray equipment now offers the best selection in flat panel digital technology.

Purchase a complete digital system or upgrade your existing x-ray system with a digital flat panel.



Digital DX-C

The Canon panel is removable and can be used for horizontal beam, weight bearing, oblique and cross table exams.



Images Direct to your Computer in seconds

Our Flat Panel digital technology delivers images in seconds, and offers lower exposure time than traditional film systems. With the Image Evaluation software, you can adjust images and diagnose immediately to inform clients in real time.

SingleTouch Monitor System

If you purchase a complete new system or upgrade from a VetRay (Sedecal) film x-ray machine, you will receive the benefit of our SingleTouch Monitor. The Single Touch Monitor is the only monitor needed in the x-ray room. Learn more at www.VetRay.com

4-way float top table

A 5 foot 4 way float top allows for easy positioning of the patient. Designed specifically for veterinarians it has a longitudinal float of 22 inches and a cross table float of 7 inches. A 6 foot table top is also available if needed.

High Frequency Generator


HF - 32 kW, 400 mA, 125 kVp, 500 mA single or three phase operation. Generator fits inside table to optimize space in x-ray room.

Hands Free Collimator Light Control

Allows operator to keep hands on patient for convenient positioning.

We are the Manufacturer

The VetRay line of x-ray solutions is created by the world's leading x-ray manufacturer, Sedecal. Our team of experts and remote support staff is always on duty to assist you.

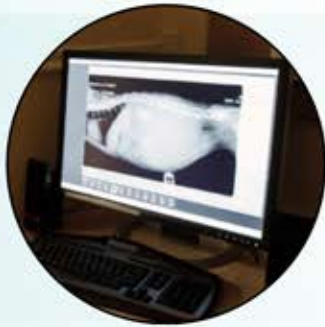

Share and transfer images on multiple platforms.

Save images in multiple formats and share with specialists immediately via DICOM transfer, burn to CD, or send to web enabled device (smartphone or tablet).



Digital DX-C

Image field size	14" x 17"
Spacial resolution	3.1 lp/mm
Scintillator	Gadolinium (Cesium also available)
Number of Bits	12 bits (14 bits A/D)
Shades of gray	4,096 (16,384)

1 year warranty if used outside of table.
3 year warranty if used stationary in table.

The Vetray system also includes:

- Acquisition PC
- 2mp high resolution 21" monitor
- Unlimited viewing licenses
- 5 year warranty on x-ray system (see detector warranty above)

800.920.9525
Info@VetRay.com



800-356-3388
978-374-6371
Fax - 978-521-2214

49 Newark Street
Haverhill, MA 01832
sales@associatedxray.com
www.associatedxray.com

www.VetRay.com

